




Double-Coated Medical Film Tape

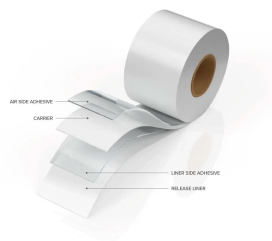

- Backing Carrier:** 0.5 mil polyester film
- Adhesive AS/LS:** Non-sensitizing acrylic (M106)
- Release Liner:** 76lb polycoated bleached kraft, silicone treated on both sides

Applications

- Surgical drape fabrication and assembly
- Medical device assembly
- Medical device attachment
- Surgical drape attachment: drape to drape, drape to skin

Features

Benefits

PET Film Carrier

- Supports for aid in processing and application
- Conformable

Non-Sensitizing Acrylic Adhesive

- High tack for quick stick application
- Suitable for use on patient's skin

76lb PC Bleached Kraft

- Low COF for hand apply applications
- Good dimensional stability and moisture resistance

Typical Test Values

	Standard UOM	Metric UOM	Test Method
Total Thickness (with Liner)	8.7 mils	0.221 mm	PSTC 133
Carrier Thickness	0.5 mils	0.013 mm	PSTC 133
Adhesive Thickness (Air Side)	1.3 mils	0.033 mm	PSTC 133
Adhesive Thickness (Liner Side)	1.3 mils	0.033 mm	PSTC 133
Liner Thickness	5.6 mils	0.142 mm	PSTC 133
Adhesion to Steel (Air Side)	38 oz/in	4.16 N/cm	PSTC-101
Adhesion to Steel (Liner Side)	48 oz/in	5.25 N/cm	PSTC 101 Mod
Recommended Storage Temperature	60 to 80 °F	15.6 to 26.7 °C	
Recommended Storage Humidity	40 to 60 %		

Values shown are based on laboratory testing and represent typical product characteristics. Actual performance will vary depending on application conditions, and no guarantee of performance is made. The user is responsible for determining suitability for the intended application. No express or implied warranty, including merchantability or fitness for a particular purpose, is made.